



Introducing the ProForm NX™

The NeXt Generation

ProForm NX

- A cost effective, customizable solution from an off the shelf product
- For clients with challenging positioning and postural support requirements
- Incorporates Air-Foam Floatation Technology
- Easy to modify components using everyday tools
- Lightweight
- Description of ProForm NX components
- Procedure for modifying ProForm NX




Contoured Base



Contoured Base

- Easily modifiable with standard, everyday tools
- Adductor Laterals
- Medial Thigh Separator
- Pre-ischial Bar
- Deep Ischial Pan
- Flex Hinge
- Modification Guides

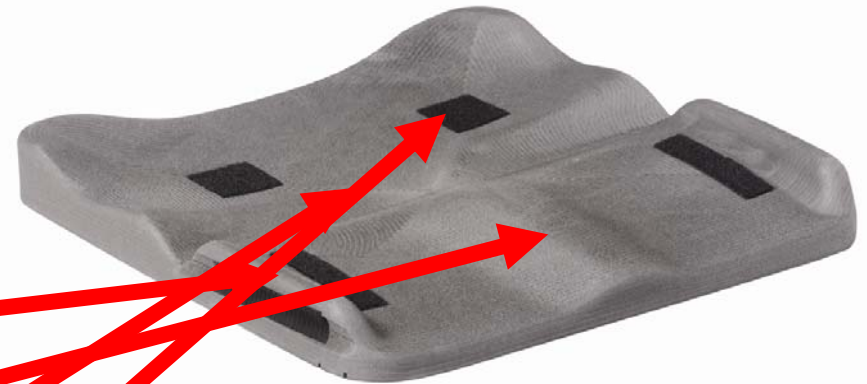



Modifiable Wedge





Modifiable Wedge

- Cut to any length for desired amount of “dump” or “squeeze”
- Because the Contoured Base is hinged, the seat-to-back angle does NOT change when Modifiable Wedge is utilized
- Can be placed inside the cushion cover





Air-Foam Floatation Cushion





Air-Foam Floatation Cushion

- Self-inflating, air and foam cushion
- Two-Way Air Valve releases air to immerse client into the cushion



Air-Foam Floatation Cushion

- Dual chambered cushion
 - Patent pending technology
 - Provides accommodation of pelvic asymmetries
- Single chambered cushion
 - Provides symmetrical pelvic and postural support





Modifiable Thigh Cushion



Modifiable Thigh Cushion

- Reticulated foam provides postural support and comfort through the thigh region
- Designed to stretch and contour over the medial thigh separator and adductor laterals



Modifiable Thigh Cushion

- Patterned for easy modification in 1" and 2" increments
- Mesh covering allows for moisture to migrate through the foam and away from skin





Cover





Cover

- Available in both mesh and incontinence
- Machine washable
- Non-skid surface on bottom



Cover

- Reticulated foam inside promotes air exchange and improves pressure distribution
- Hook and loop aligned on the bottom allows for cinching of the cover when making modifications



Modifying the ProForm NX

- Modify Contoured Base and Modifiable Wedge using common blade or electric knife
- Cinch cover using provided hook and loop
- Alter Modifiable Thigh Cushion with scissors





ProForm NX is the seating solution for:

Pelvic Obliquities





ProForm NX is the seating solution for:

Upper Leg Length Discrepancies





ProForm NX is the seating solution for:

Wind Swept Deformities





ProForm NX is the seating solution for:

CVA/ Unilateral Scooters





ProForm NX is the seating solution for:

Above Knee Amputees (AKA)





ProForm NX is the seating solution for:

Clients prone to sacral sitting





ProForm NX is the seating solution for:

Symmetrical sitters requiring extra support





ProForm NX is the seating solution for:

Clients with limited hip flexion





ProForm NX is the seating solution for:

Clients with hot spots



Conclusion

- A cost effective, customizable solution from an off the shelf product
- For clients with challenging positioning and postural support requirements
- Approved for CMS Code K0736 Skin Protection and Positioning
- Available in sizes 12 x 14 to 20 x 20

